

BELLIGER ACE



ALL-IN-ONE
MOBILE FPD C-ARM


A New Standard for C-ARM: Integrated FPD



All-in-one design



3.2kW HFG

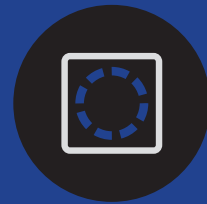

201×200mm CMOS /
21×21mm a-Si TFT FPD



Extended
free space of 860mm



15.6" Touchscreen panel



Digital Subtraction
Angiography




201×200mm CMOS / 210×210 mm a-Si TFT-based FPD, with acquired images viewable on a 4K UHD / FHD 27" monitor.

All-in-one monitor design enables easy installation and mobility, while providing a spacious 860mm procedural area, enhancing convenience for both medical staff and patients.

**Guaranteed Image Quality
from a Specialist in C-Arms**

Clinical Data

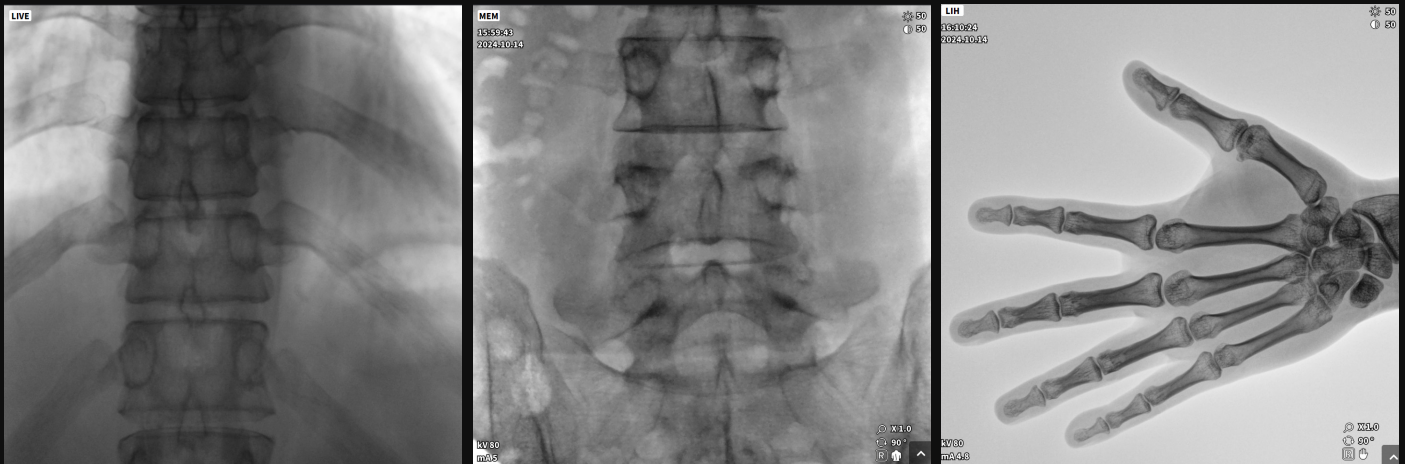

Pain Management case



Live Image with contrast medium injected

Dimension

Depth : 1,865 mm


Width : 800 mm

Height


- At moving Monitor Arm :

1,808(Stadard) - 2,168(Max.) mm

Weight : 310kg

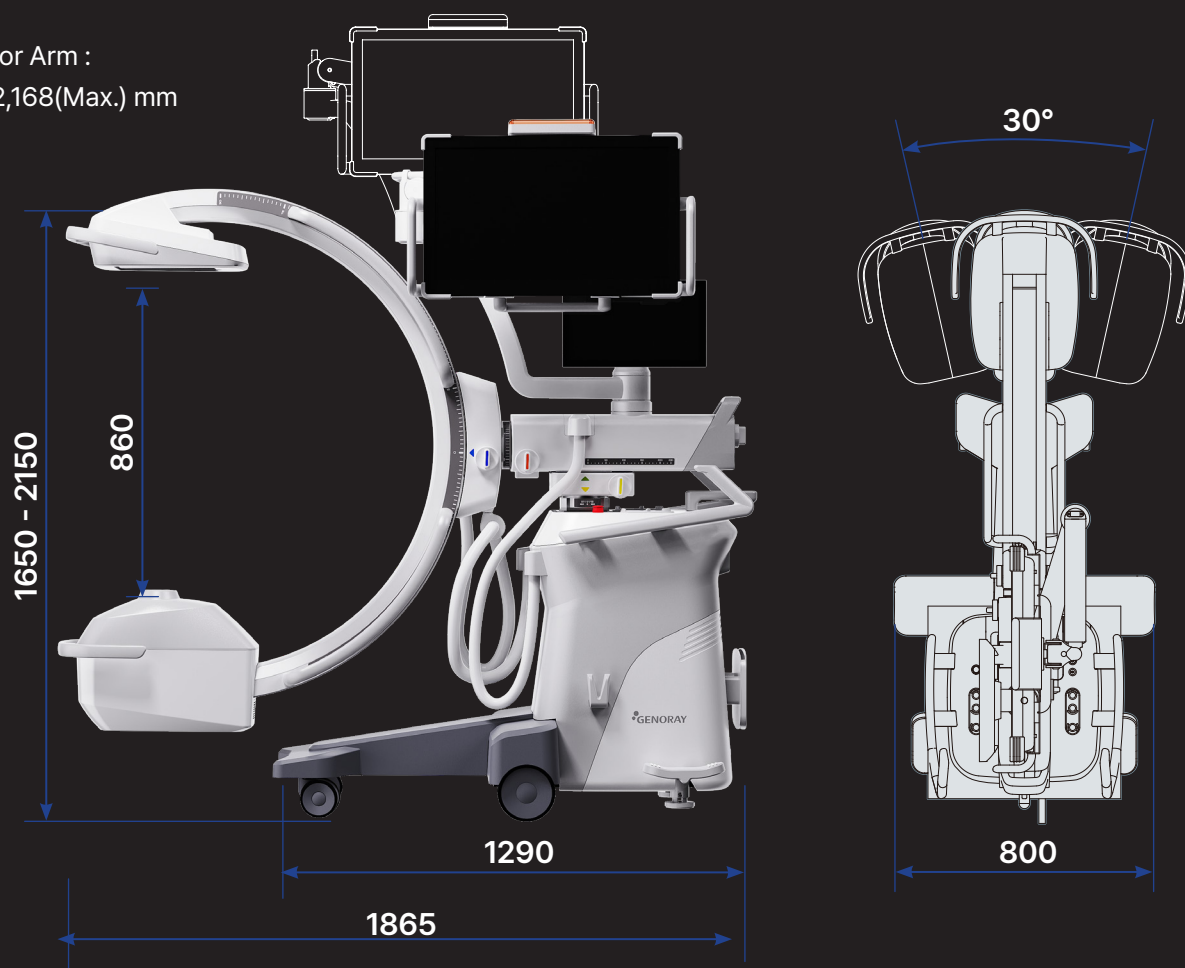

Room size

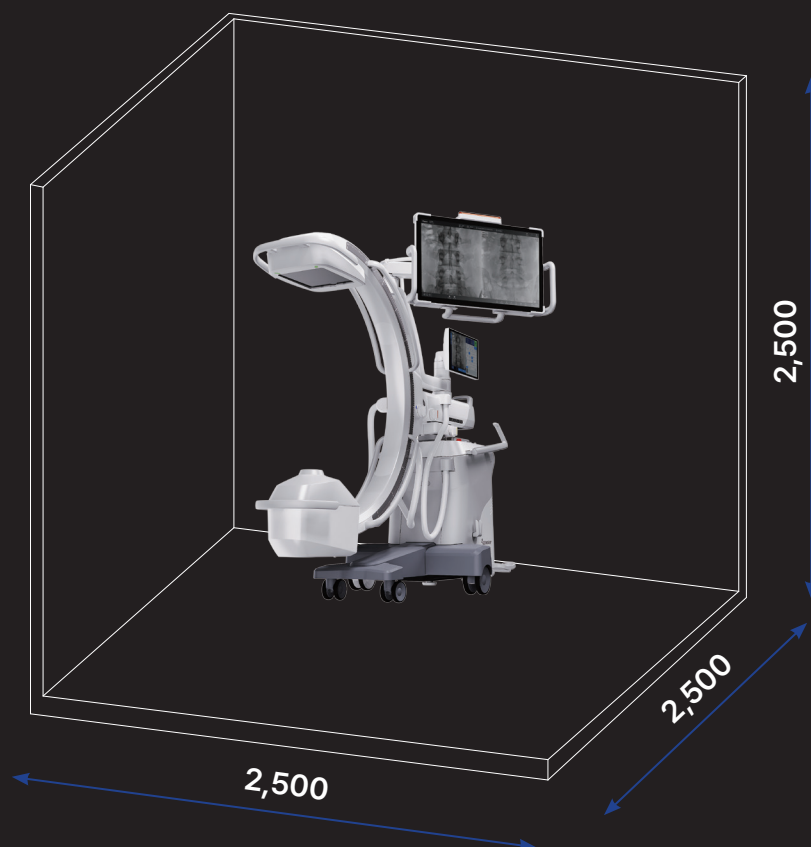

Support Your Expertise,

ZENIS IV plus





Enhanced software

Proprietary C-Arm Software
Developed by GENORAY



15.6" Touchscreen Interface

A touch interface that allows easy control of system operations and image processing.



ABC (Auto Brightness Control)

Automatically searches for the appropriate exposure conditions for the X-ray target area, delivering optimized real-time imaging.



Motion Correction Function

Minimizes afterimages caused by patient movement during image acquisition.



Metal Correction Function

Adjusts metal artifacts during surgical tool imaging to provide optimal images.



Style Template Function

View acquired images in three versions using pre-configured style templates.

Experience Advanced Imaging Technology with **Enhanced Software Features**

ZENIS IV plus



GENORAY Co.,Ltd.

T +82-31-5178-5500
E genoray@genoray.com
W www.genoray.com

GENORAY EU GmbH

T +49-30-50969498
E info@genorayeurope.com
W www.europe.genoray.com

GENORAY Türkiye

T +90-216-463-57-37, +90-216-463-57-53
E genorayturkiye@genoray.com
W www.turkiye.genoray.com

GENORAY America Inc.

T +1-855-436-6729
E inquiry@genorayamerica.com
W www.genorayamerica.com

GENORAY Japan

T +81-45-620-4971
E Info@genorayjapan.co.jp
W www.genorayjapan.co.jp

GENORAY Shanghai

T +86-021-6785-8538
E genoraychina@genoray.com
W www.genoray.com.cn